



Innovative Technologie Völp



TBI STEP

PRODUCT INFORMATION

www.it-v.net


WHY TBI STEP

The homogeneous dose distribution with simultaneous maximum protection of healthy tissue and organs at risk represents one of the greatest challenges in total body irradiation (TBI) in VMAT- technology.

From a height of approx. 120 cm, this is not possible in one patient orientation (e.g. HFS) for technical reasons. This means that a mandatory change to the FFS is necessary.

The **TBI STEP** enables this change without the active help of the patient by a rotating the table top by 180° with the help of the appropriate adapter and handles on the Linac table. This makes it possible to carry out the irradiation in an acceptable period of time with the greatest possible patient comfort with a maximum of 5 isocentres and without the risk of field overlaps.

An individual, external shielding of the OARs is unnecessary due to the VMAT- technology. The integrated indexing system offers optimally reproducible positioning of the patient by means of common positioning devices.

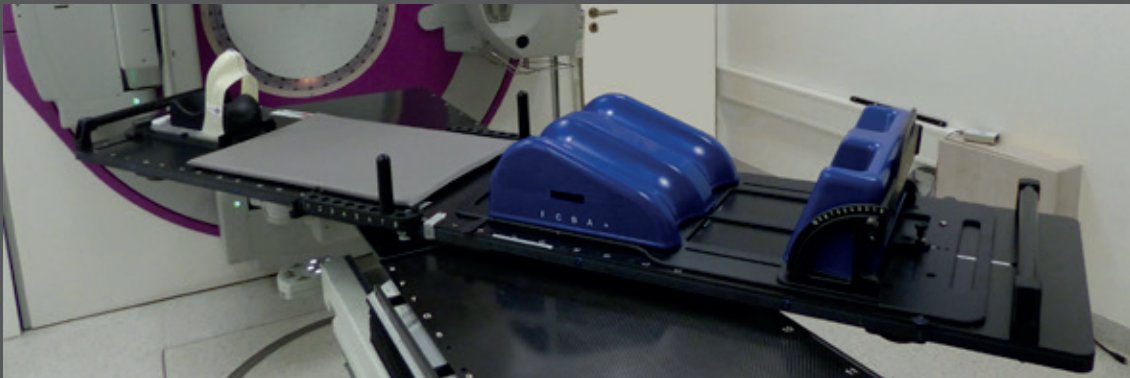

TBI STEP

Rotatable Table-Top for TBI in VMAT



- » homogeneous dose distribution over the whole PTV by using VMAT-technology
- » best possible protection of healthy tissue and OARs without using external shields
- » no active change of the patient from HFS to FFS necessary
- » simple rotation of the **TBI STEP** in a few seconds by means of handgrips
- » Minimization of unexpected patient movement through use conventional positioning devices
- » precisely reproducible positioning of the patient
- » Avoidance of incorrect isocenter shifts
- » Avoidance of gantry-couch collision by using permitted Isocentric coordinates

System TBI STEP Komponenten:

P10107-615 CT STEP Basic Module 2,40m

P10107-673 CT STEP Basic Module 2,40m for Varian

P10107-527 TBI Adapter CS

P10107-409 Indexing Bar 14 Carbon (set of 2)

P10107-511 Hand Grip Pole (left & right)

IT-V Innovative Technologie Völp

Exlgasse 20a
6020 Innsbruck, Österreich
T +43 664 411 33 24
F +43 512 27 75 77
office@it-v.net
www.it-v.net

Service & Support International | Back-Office

Elisabeth Entner
T +43 664 411 33 24
e.entner@it-v.net

Service & Support Germany | France | CH

Boris Jäger
T +49 151 275 791 99
b.jaeger@it-v.net

Service & Support Germany | NL | UK

Kerstin Grings
T +49 172 29 83 887
k.grings@it-v.net

Service & Support Germany North-East

Rudi Schröder
T +49 171 55 44 690
r.schroeder@it-v.net

Service & Support Germany | CZ | SK

Anja Munk
T +49 172 617 85 12
a.munk@it-v.net

www.it-v.net